



Features:


- Medical-grade HD CMOS sensor to ensure color saturation and image authenticity
- Ergonomic design of double-side freeze button
- seamless welded shell, easy to clean and disinfect
- High luminance LED make the images much clearer

# SuperVista<sup>TM</sup> HD

Auto-focus (fixed) Intraoral Camera HDI-210

Working Scheme :



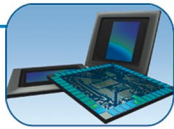

#### Industrial Image sensor

1.3 megapixel HD CMOS sensor could ensure high quality intraoral images.



#### Professional dental lens

Professional lens is an important precondition to generate high quality images which reflect the objects truly and make it easier to take photos.



#### Direct digital imaging

USB2.0 interface with direct digital imaging, no capture card required, achieve higher speed and lossless images.



#### Powerful imaging management software

HandyDentist provide multifunction for image processing and case management which facilitate better communication between dentists and patients.



#### Support the third party software

SDK of hardware are available to OEM customers and our Twain driver make the hardware possible to work smoothly in any other software with Twain function.



#### Double-side freeze button

Ergonomic design of double-side freeze button brings dentists much more convenience in practice

### Specification

Image Sensor	1/4" HD CMOS
Active Pixels	1280 x 1024 (1.3MP)
Resolution	1280 x 720 ( 720p )
Frame Rate	30fps@720p
Focus Range	5mm~35mm
Angle of View	80°
Distortion	<5%
Lighting	6 LEDs
Signal Output	USB 2.0
Length of Cable	2.5m
Driver Type	UVC/WDM
TWAIN	Yes
Operation system	Windows 2000/XP/Vista/Windows 7/Win8(32bit&64bit)

## Shanghai Handy Medical Equipment Co.,Ltd

Tel : +86 021-56766108 Fax : +86 021-56766308

Website : [www.handycreate.com](http://www.handycreate.com)

E-mail : [h.sale@hotmail.com](mailto:h.sale@hotmail.com) [handy@handycreate.com](mailto:handy@handycreate.com)

After - sales supports:

Skype: [handysupports](https://www.skype.com/handysupports) Email: [support@handycreate.com](mailto:support@handycreate.com)